


HOME CARE WITH THE GRASBY SYRINGE DRIVER

You are going home with a Grasby Syringe Driver to pump your medication dosage through as ordered by your doctor.

Here are a few tips to help you care for it, and to help when troubleshooting problems.

Measuring
Scale

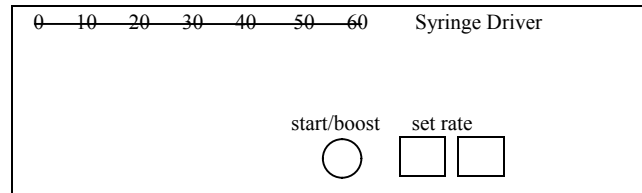

Start Button

Setting

ALARM	PUMP ACTION	REASON
Continual beeping	Stops	Syringe is empty - requires new syringe of medication
Continual beeping	Stops	Tubing is blocked - may be kinked - requires unkinking
If other reason - may require new needle for access - ring the hospital and come in.		
Light stops flashing	Low battery - needs replacing.	
You may replace it yourself if you have a spare - after replacement you will need to press the start/boost button - check for flashing light. You may choose to return to hospital for a replacement battery.		

HINTS: Do not get your syringe driver wet when bathing.

Keep your driver in the case, and in the given pouch or a pocket of your clothing.

Check the pump regularly (ie. Every 1 - 2 hours to ensure the light is flashing - this signals that the pump is working.

If you are having any problems, or are unsure of what to do next, please ring the nursing staff - Day Oncology, Mon - Fri, 8.00am to 5.00pm - 3840 195

Ward 10B - at other times - 3840 8973.